 PT ETANA BIOTECHNOLOGIES INDONESIA	STANDARD OPERATING PROCEDURE		Page : 1 of 8	
	TITLE : PERAWATAN TERHADAP MESIN PARTICLE COUNTER MAINTENANCE OF PARTICLE COUNTER MACHINE		SOP No. : SOP/EBI/EN-063 Revision : 00 Effective Date : Reviewed Date :	
PREPARED BY :	REVIEWED BY :	APPROVED BY :	APPROVED BY :	
WENDI RUKMANSYAH (ENGINEERING SPV)	PURNO BUDI KISWANTO (ENGINEERING MANAGER)	VINO SOADUON (PROD MANAGER)	SHANDY (QA MANAGER)	

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Provide clear instructions on how to maintenance for Particle Counter Machine in the production area.

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This SOP is used by the Engineering department as a guideline for carrying out Particle Counter Machine maintenance in the production area.

3. TANGGUNG JAWAB

RESPONSIBILITY

3.1. Teknisi Engineering *Preventive Maintenance* bertanggung jawab :

Engineering Technician Preventive Maintenance are responsible :

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Performing maintenance duties properly, in accordance with this SOP.

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Perform maintenance according to a predefined schedule.

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Report any damage and repairs that are needed to the Engineering Supervisor.

3.2. Supervisor Engineering bertanggung jawab :

Engineering Supervisor are responsible :

3.2.1. Membuat jadwal *maintenance* secara berkala.

Make schedule of maintenance periodically.

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Ensure implementation of maintenance is according to the schedule.

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Ensure all measuring instruments involved in the maintenance process are in accordance with the specifications specified (calibrated).

3.2.4. Mengadakan *review* bulanan terhadap hasil kerja Teknisi *Preventive Maintenance*.

Perform a monthly review of the workings of the Preventive Maintenance Technicians.

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Provide training to others department which related.

3.3. Manajer Engineering bertanggung jawab :

Engineering Manager is responsible :

3.3.1. Menentukan *spare parts* pengganti apabila *spare parts* yang sesuai sudah tidak dijual lagi.

Determine replace parts if suitable spare part are no longer sold.

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Resolve and coordinate necessary actions, if necessary replacement, equipment modification, system changes and prepared change control.

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The Production Machine Operator is responsible for :

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Report if the area is detected by defects / abnormalities that normally impede the operation of the machine.

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Perform daily maintenance such as: tighten ferulle clamps loose bolts, clean the body of the machine / Main Frame at all times so that it does not rust and the surface of the machine remains good, cleaned every machine that is the responsibility, especially on the inside of the crust that is attached.

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Ensure the condition of electrical and electronic equipment that is part of the machine placed in water-prone areas, in order to always avoid splashing water.

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Conduct condensate water discharge at Air Service Regulator on machine equipped with pneumatic system.

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The Calibration Staff of the Quality Assurance Department is responsible for calibrating the components of the measuring instruments used, such as Thermometer, Pressure Gauge, Thermo Controller etc., and conducting regular qualifications..

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The particle counter machine works by detecting particles that pass through the sensor and counting them based on the particle size. This helps in ensuring that a particular environment or product meets the set hygiene standards.

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Autonomous maintenance is an independent engine maintenance carried out by the machine operator. If all this time the operators are only trained to operate the engine, then it's time to be further trained.

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Before implementing preventive maintenance, engineering staff must coordinate with the operator about the area to prepare the area.

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When carrying out maintenance, the related operator must accompany Engineering to avoid operating errors when the machine is running or the operator can undergo autonomous maintenance on the machine.

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Work safety must be paid attention to at the time of preventive maintenance, use PPE (Personal protective equipment) which according to allotment in condition when maintenance.

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For the replacement of critical components or spare parts that have a direct impact on the product when replaced with different brands and different types but one type, make change control.

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1 set of toolbox, when performing Preventive Maintenance.

6.2. BAHAN

MATERIAL

6.2.1. Ethanol 20%

Ethanol 20%

6.2.2. Kain Wipes

Wipes Fabric

7. PROSEDUR

PROCEDURE


Perawatan 3 bulanan (L2) dan 1 Tahunan (L4) pada Mesin *Particle Counter*
Maintenance 3 month (L2) and 1 Year (L4) of Mesin Particle Counter Machine.

7.1. Perawatan 3 bulanan (L2).

Maintenance 3 month (L2).

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Maintenance of Electrical elements.

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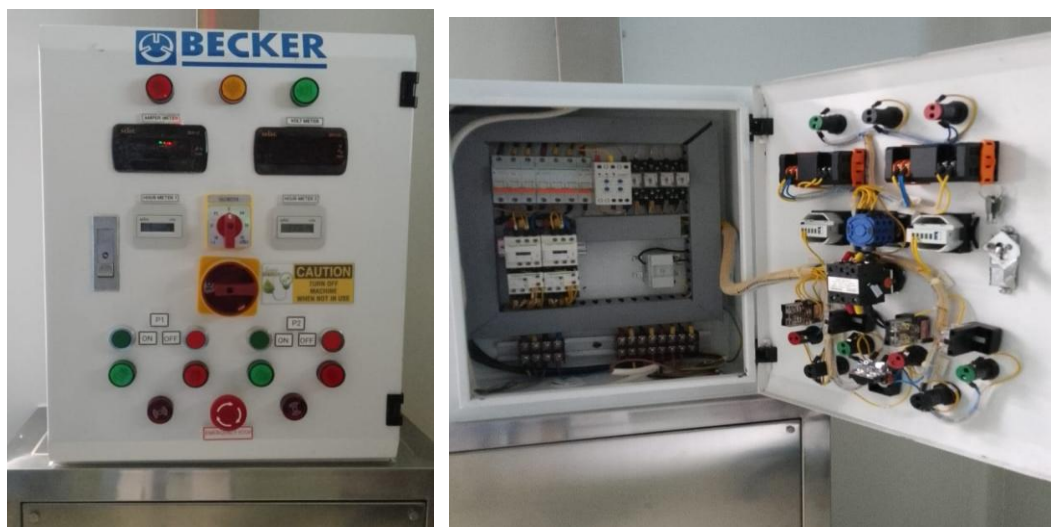Figure 1. Particle counter machine electrical panel for bosch filling machine.

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Check electrical power connections 3 Phase and other connections.

7.1.1.2. Periksa tombol - tombol pengoperasian pada panel.

Check the operating buttons on the panel.

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Clean the panel from dust using a brush.

7.1.1.4. Periksa kondisi MCB, jika terlihat ada yg rusak ganti MCB nya.

Check the condition of the MCB, if you see anything damaged, replace the MCB.

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Check and make sure the bolts securing the skun are tightened.

7.1.2. Perawatan pada unit *vacuum pump becker*.

Maintenance on the Becker vacuum pump unit


Gambar 2. Unit mesin vacuum untuk particle counter mesin filling bosch

Figure 2. Vacuum machine unit for Bosch filling machine particle counter

- 7.1.2.1. Bersihkan *body unit vacuum pump* yang berada di ruangan *vacuum pump*.
Cleaning body unit vacuum in the vacuum pump room.
- 7.1.2.2. Bersihkan filter udara pada unit *vacuum*.
Cleaning air filter in the unit vacuum.
- 7.1.2.3. Periksa kondisi *valve – valve* pada *header vacuum*, ganti jika *valve* terlihat rusak atau bocor.
Check the condition of the valves on the vacuum header, replace them if the valves appear damaged or leaking.
- 7.1.2.4. Periksa sambungan – sambungan pada selang mulai dari unit *vacuum* menuju ke ruangan mesin *filling bosch*.
Check the connections on the hose starting from the vacuum unit to the Bosch filling machine room.
- 7.1.2.5. Periksa sambungan pada *probe vacuum* yang berada di bagian dalam mesin *filling bosch*, kencangkan koneksinya jika terdapat yang kendur.
Check the connections on the vacuum probe inside the Bosch filling machine, tighten the connections if there are any loose ones.
- 7.1.2.6. Periksa tekanan *vacuum* yang terbaca pada *pressure gauge*, tekanan *vacuum* normal adalah – 0,70 sampai -0,80Mbar.
Check the vacuum pressure read on the pressure gauge, normal vacuum pressure is – 0.70 to -0.80 Mbar.

7.2. Perawatan 1 Tahunan (L4).

Maintenance 1 year (L2).

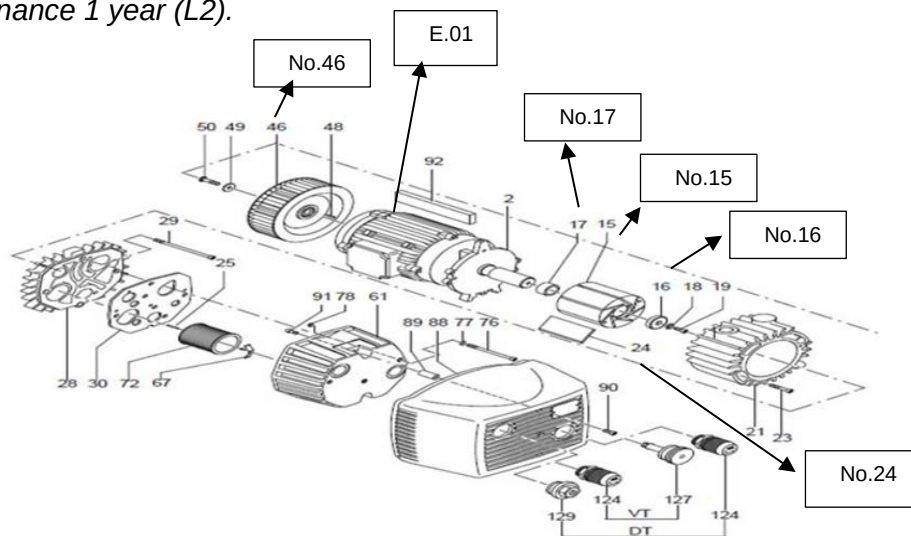Gambar 3. Vacuum Pump Becker

Figure 3. Becker vacuum pump

7.2.1. Perawatan pada unit vacuum pump becker.

Maintenance on the Becker vacuum pump unit.

- 7.2.1.1. Periksa *bearing vacuum*, ganti *bearing* jika putaran vakum terdengar kasar diputar secara manual sedikit kasar, contoh bearing terdapat di nomer 16,17 dan bagian motor listriknya Lihat pada gambar 3. Vacuum Pump Becker. Point E.01

Check the vacuum bearing, replace the bearing if the vacuum rotation is roughly rotated manually a little rough, an example of the bearing is found at number 16, 17 and the electric motor part See Figure 3. Becker Vacuum Pump. Point E.01

- 7.2.1.2. Periksa kipas pendinginan pada *vacuum*, kencangkan jika kendur, contoh fan pendingin terdapat di nomer 46. Lihat pada gambar 3. Vacuum Pump Becker. *Check the cooling fan on the vacuum, tighten if loose, an example of a cooling fan is found at number 46. See Figure 3. Becker Vacuum Pump.*

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Check the connection of the electric motor to the carbon blade housing tighten it if it is loose or replace it if the condition looks bad, an example of the connection of the electric motor to the carbon blade housing is number 15.

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Check the condition of the carbon blade, replace it if the condition of the carbon blade looks damaged or the size changes from the standard size, the standard size adjusts to the carbon blade housing distance, an example of a carbon blade is found at number 24. See Figure 3. Becker Vacuum Pump.

8. REFERENSI

REFERENCE

Buku manual perawatan berkala Mesin Penghitung Partikel
Manual book of periodic maintenance Mesin Particle Counter

9. LAMPIRAN

ATTACHMENT

Formulir perawatan berkala terhadap Mesin *Particle Counter* SOP/EBI/EN-063-F01A.
Periodically Maintenance form of Particle Counter Machine List No. SOP/EBI/EN-063-F01A.

10. RIWAYAT

HISTORY



Tanggal <i>Date</i>	No. Change Control <i>Change Control No.</i>	No. Revisi <i>No. Revision</i>	Alasan Perubahan <i>Reason of Change</i>
	CC-EBI-EN-018-24	00	Dokumen baru <i>New Document</i>

11. DISTRIBUSI

DISTRIBUTION

- 11.1. Departemen Engineering.
Engineering Department.
- 11.2. Departemen Produksi.
Production Departement
- 11.3. Departemen Sistem Mutu
Quality System Departement
-

 PT ETANA BIOTECHNOLOGIES INDONESIA	FORMULIR / FORM								Document No. : SOP/EBI/EN-063-F01A	
	PERAWATAN TERHADAP MESIN FUME HOOD MAINTENANCE OF FUME HOOD MACHINE									

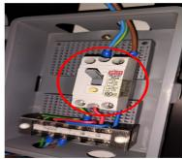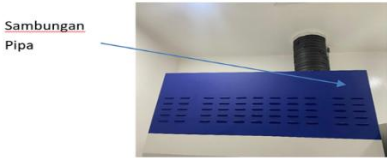Level :											
Month :											
Date :											
Day :											
Hours meter :											
Level	Component	Activity							Condition		Remark
		C	I	L	R	S	T	W	B	G	
L2											
	Rotation of the blower motor										
											
	Ampere on the blower motor										<2A
											
	Connection the blower pipe										
											
	1 phase electrical connection										
											
	Electrical installation										
											
	valve connection on the soft water supply										
											
	Hinges on the door										

 PT ETANA BIOTECHNOLOGIES INDONESIA	FORMULIR / FORM		Document No. : SOP/EBI/EN-063-F01A
	PERAWATAN TERHADAP MESIN FUME HOOD MAINTENANCE OF FUME HOOD MACHINE		

Note :

No	Date	Start	Finish	Duration		Total personnel	Area	Remark
				Hour	Minute			

Remark : C : Cleaning (pembersihan) I : Inspection (pemeriksaan) L : Lubrication (pelumasan) R : Replacement (penggantian) S : Setting (penyetelan) T : Tightened (dikencangkan) W : Welding (pengelasan) G : Good (baik) B : Bad (jelek)	Level : L2 : Every 3 month	Done by : ()	Checked by : ()
--	-----------------------------------	---	--